

Simplify your workflow
without sacrificing bond strength.



Nexus™ Universal

A Truly Universal Adhesive Resin Cement

Nexus Universal eliminates the need for primers or activators – cutting the complexity of typical cementation processes by up to 50%. With simplicity comes added flexibility and convenience. Nexus Universal adapts to your preferred cementation technique, so you don't have to change your protocol. It also pairs with leading universal bonding agents including OptiBond™ Universal, part of the OptiBond™ family of products widely recognized as the gold standard.

Features & benefits



Simplify Cementation

Nexus Universal reduces the complexity of cementation steps by up to 50%. You'll enjoy a simplified, and more efficient workflow with fewer opportunities for error.



Do it Your Way

You'll have the freedom to keep using your preferred cementation protocol and stay confident in a strong bond – with reduced steps.



Better Bond Strength

Nexus Universal delivers one of the most secure holds compared to other leading adhesive resin cements, regardless of which leading universal bonding agent you use. See charts for universal bonding agents tested.*




Color Cleanup Indicator with One-Peel™ Cleanup

Nexus Universal Chroma dispenses pink and fades to a cement shade, allowing you to recognize the precise moment your excess cement is ready to be removed. Then, with just One-Peel™, your restoration is secure.

Simplified workflow

Nexus Universal cuts the complexity of typical cementation processes by up to 50%. Simple Dark Cure Technique with any OptiBond adhesive or any other leading universal bonding agent.

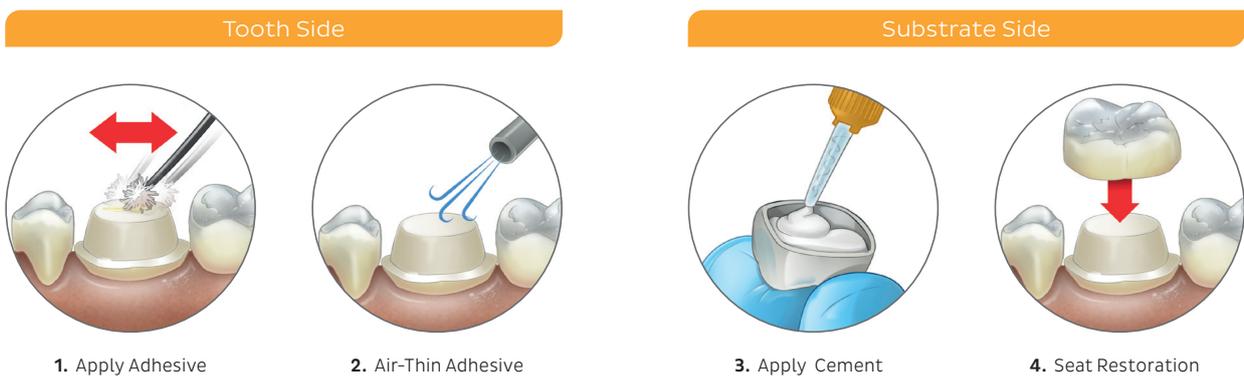


- No need for a dual-cure activator
- No need to light cure the adhesive
- No need for metal primers / silane



A simplified, faster, and more efficient workflow means less opportunities for error.


Boost your bond strength

Nexus Universal and OptiBond Universal deliver stronger bonds than other leading cementation systems in dark-cure mode.

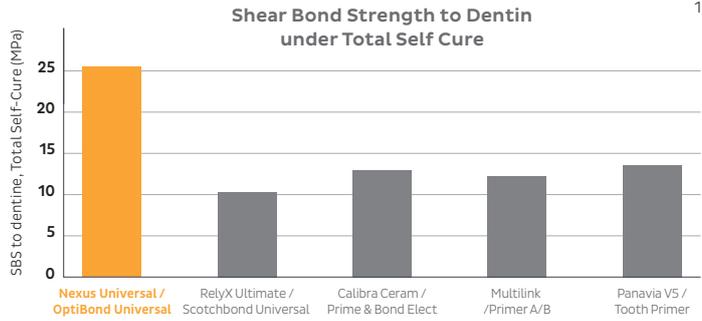


Experience more color stability

Nexus Universal features a proprietary redox initiator system which virtually eliminates discoloration common in BPO/tertiary initiator systems. This results in better color stability and long-term esthetics.

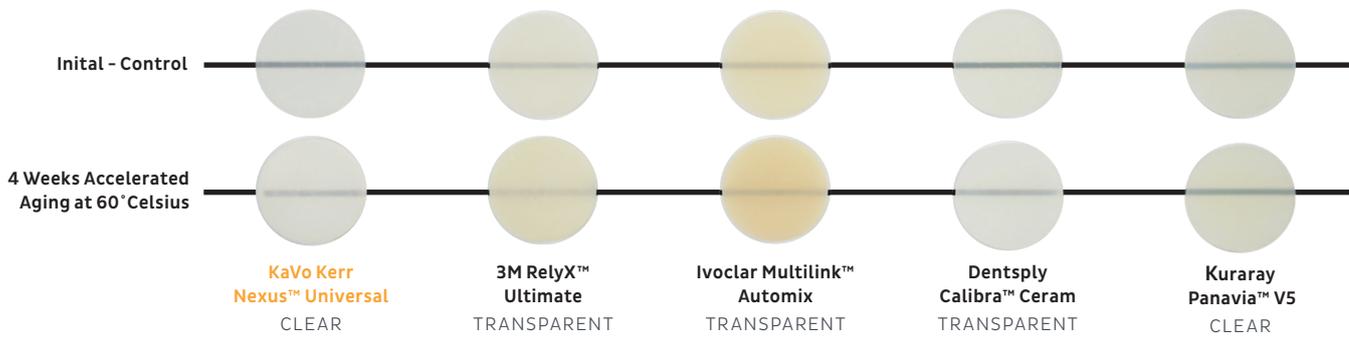


Choose your bonding agent

Nexus Universal offers total dark-cure compatibility with any leading universal bonding agent including OptiBond Universal and other OptiBond adhesives. Be confident in your bond strength, regardless of which bonding agent you use. No dual-cure activators are needed for dark-cure mode.

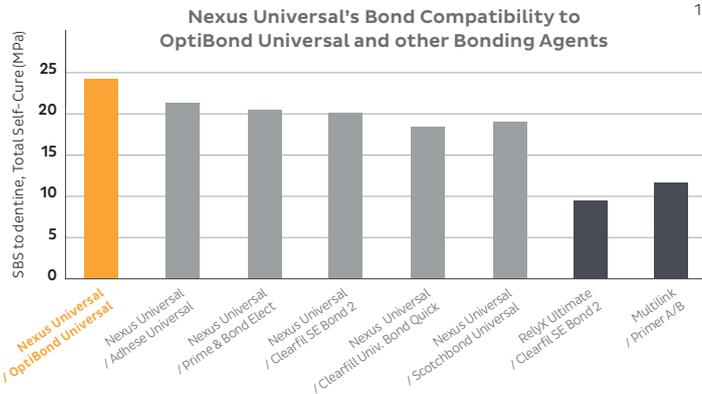

High performance on any substrate

Nexus Universal is compatible with various indirect restorations and CAD/CAM materials, including resin, ceramic, zirconia, gold, base metal alloys, alumina, lithium disilicate (e.max) and titanium.

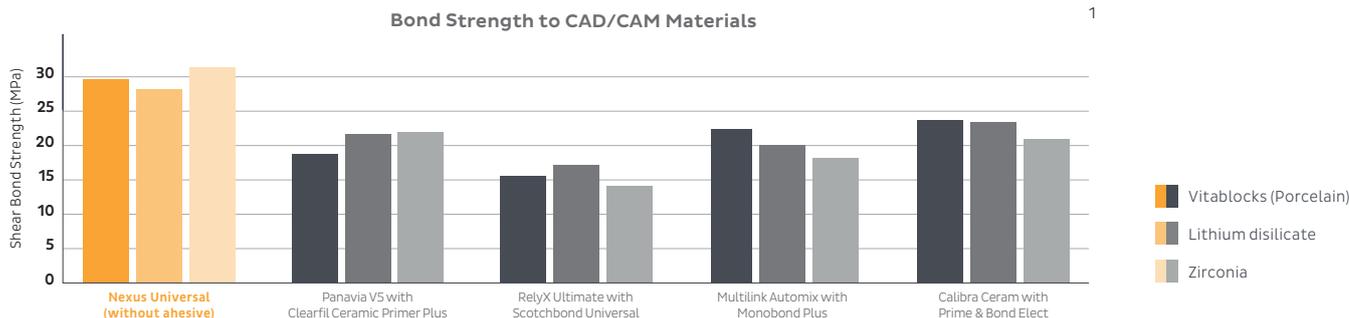


1 - Internal data available upon request

Nexus Universal Chroma color cleanup indicator

When the **pink** disappears, you're ready to cleanup!



Cement is pink when dispensed and gradually fades to the cement shade. Tack-cure or self-cure for more control over single or multiple crowns. Cleanup is simple once the pink color has turned to the cement color, indicating the optimal window of time to cleanup excess cement - reducing the risk of leaving excess cement.

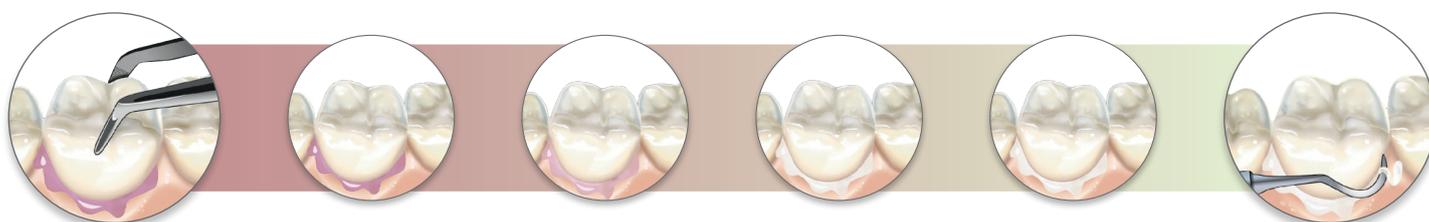

Application

Cleanup

Nexus™, OptiBond™ and One-Peel™ are either registered trademarks or trademarks of Kerr Corporation in the United States and/or other countries. All other trademarks are the property of their respective owners.

Nexus™ Universal. Refills and accessories Ordering information.

REF/Part. No.	Description
36990	Trial Kit - 1x 5g Syringe (Clear), 8x OptiBond Universal Unit Dose Pawns, Accessories
36989	Intro Kit - 4x 5g Syringes (Clear, Clear w/ Chroma, White, Yellow), 1x OptiBond Universal Bottle, Accessories
36991	Refill Kit - 2x 5g Syringes, Clear
36992	Refill Kit - 2x 5g Syringes, Clear with Chroma
36993	Refill Kit - 2x 5g Syringes, White
36994	Refill Kit - 2x 5g Syringes, White Opaque
36995	Refill Kit - 2x 5g Syringes, Yellow
36996	Refill Kit - 2x 5g Syringes, Bleach
36999	Mixing Tips - 50pk Regular
37000	Mixing Tips - 50pk Wide



24 Months Shelf Life

Easy Storage with no refrigeration required.

