




## ART<sup>3</sup>: A NEW WAY TO BULK FILL

Achieve lifelike restorations with the next generation of clinically proven Adaptive Response Technology (ART) found in Harmonize™ and OptiShade Universal Composite. Powered by Enhanced Adaptive Response Technology, SimpliShade Bulk Fill offers easy handling and low shrinkage.

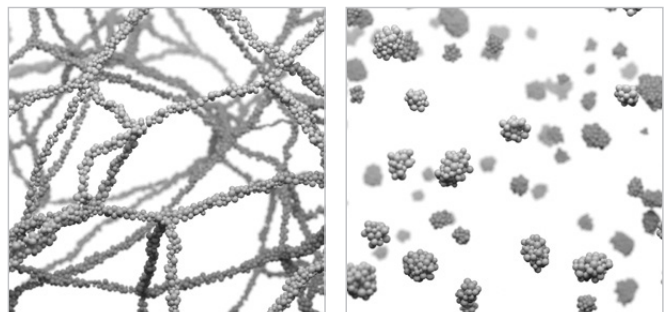


### Adaptive Rheological Response Technology for Ease of Sculpting

ART<sup>3</sup> uses rheological modifiers that reduce the composite's viscosity as you manipulate it. When you're done shaping, viscosity recovers and anatomy stabilizes before light-curing — reducing the need for resculpting.

### Adaptive Color-Blending Response Technology for Beautiful Blending

Achieve exceptional blending with special nanofillers that mimic enamel structures and interact with light in a natural way. The restoration and tooth become virtually indistinguishable.

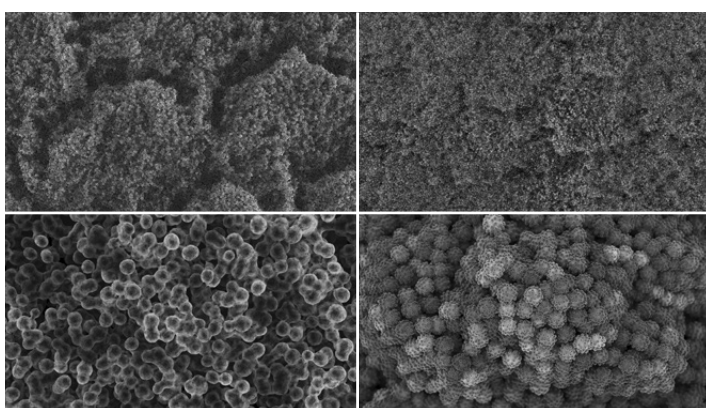


Human Enamel

SimpliShade Bulk Fill

### Adaptive Shrinkage Stress Response Technology for Durable Results

During polymerization, SimpliShade Bulk Fill uses low shrinkage resin and innovative nanofiller design, resin-infused nanofillers boost elasticity for stress relief and a special coating eases interfacial stress. Uniting three shrinkage reduction mechanisms, this enhances marginal adaptation and integrity by minimizing composite restoration stress on the adhesive layer for enhanced marginal adaptation.

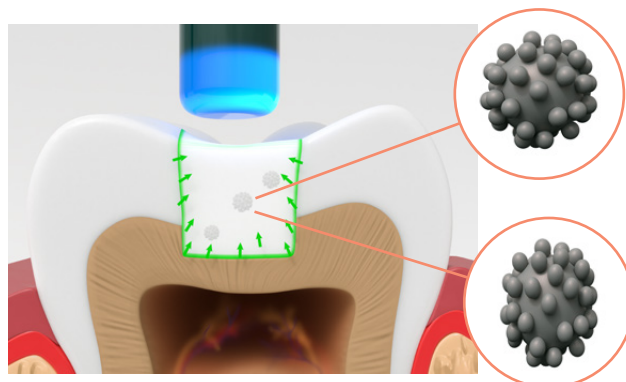

Learn more about ART<sup>3</sup> and the effect on shrinkage stress

## Ordering Information



### SimpliShade Bulk Fill

Part #	Item Description
37211	SimpliShade Bulk Fill Composite Syringe
37212	SimpliShade Bulk Fill Composite Unidose 20 Pack

## KEEP IT IN THE FAMILY

with our full line of SimpliShade products!



### SimpliShade Bulk Fill Flow

Part #	Item Description
37217	SimpliShade Bulk Fill Flow Syringe
37218	SimpliShade Bulk Fill Flow Unidose 20 Pack



### SimpliShade Universal Composite

#### Unidose

Part #	Item Description
37006	Light - 20 Pack
37007	Medium - 20 Pack
37008	Dark - 10 Pack
37009	Bleach White - 10 Pack
37010	Universal Opaque - 10 Pack

#### Syringe

Part #	Item Description
37001	Light
37002	Medium
37003	Dark
37004	Bleach White
37005	Universal Opaque



For an ultra high shine finish, use Kerr's line of Polishers!



One Shade. One Step. A Better Way to Bulk Fill.

For more information scan the QR Code or visit [go.kerrdental.com/simplishadebulkfill](http://go.kerrdental.com/simplishadebulkfill)

# SimpliShade™ BULK FILL



## Simplify with One Shade and One Step

Beautiful restorations just got simpler, with SimpliShade™ Bulk Fill. You can skip the capping layer with a one-step placement and reduce inventory requirements with one shade for all 16 VITA® classical shades. Powered by Ternary Adaptive Response Technology (ART<sup>3</sup>), SimpliShade Bulk Fill is easy to handle with low shrinkage stress.

[kerrdental.com](http://kerrdental.com)

©2023 Kerr Corporation. SimpliShade and Harmonize are trademarks of Kerr Corporation. All other trademarks are properties of their respective owners. Data on File. MKT-22-0879\_UK Rev 0

**Kerr™**

**Kerr™**



# ONE SHADE IS ALL IT TAKES

Get the adaptability, durability, depth of cure, and simplicity you need with SimpliShade Bulk Fill.

SimpliShade™ Bulk Fill is indicated for Class I-V restorations - direct placement of anterior and posterior restorations, core build-ups, and repair of defects in restorations.



## One-Shade System

Using one shade that blends seamlessly with all 16 VITA® classical shades, SimpliShade Bulk Fill saves you time selecting shades and managing inventory.



## One Step

Eliminate the need for a capping layer, complete restorations with a SingleFill™ composite. With its enhanced curing mechanism, SimpliShade Bulk Fill makes it easy to achieve polymerization depths up to 5mm in a single step.



## Simplicity with Durability

SimpliShade Bulk Fill offers lower wear rate, higher flexural strength, and comparable shrinkage stress relative to other bulk fill brands on the market, keeping your restorations stronger and long-lasting.




## Chameleon Effect


Special nanofiller architecture mimics the natural enamel structure and interacts with light for virtually indistinguishable blending.

## Flexural Strength

SimpliShade Bulk Fill demonstrates comparable or higher flexural strength compared to Competitor 1 and Competitor 2 resulting in better resistance to chipping and fracture over time.

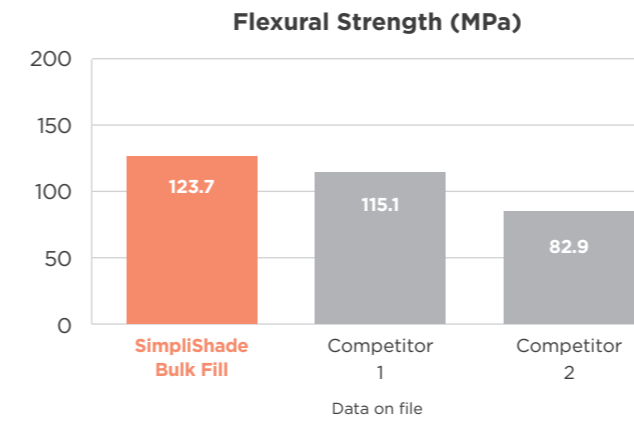

## Shrinkage Stress (at 2 min.)

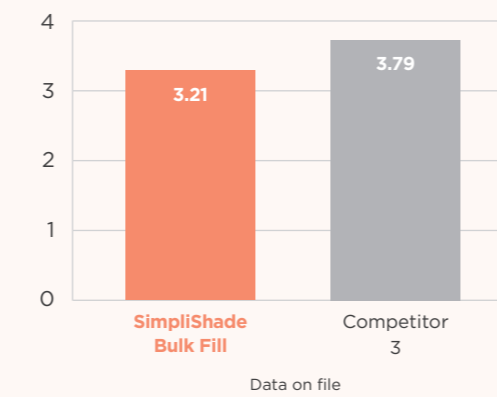


## Shrinkage Stress


SimpliShade Bulk Fill has comparable shrinkage stress to Competitor 3, resulting in lower marginal degradation and postoperative sensitivity for a long-lasting restoration.

## Wear Rate

SimpliShade Bulk Fill has the lowest wear rate relative to Competitor 3, Competitor 1, and Competitor 2. Low wear rate results in a restoration that will withstand the daily forces of chewing, brushing, and other occurrences.

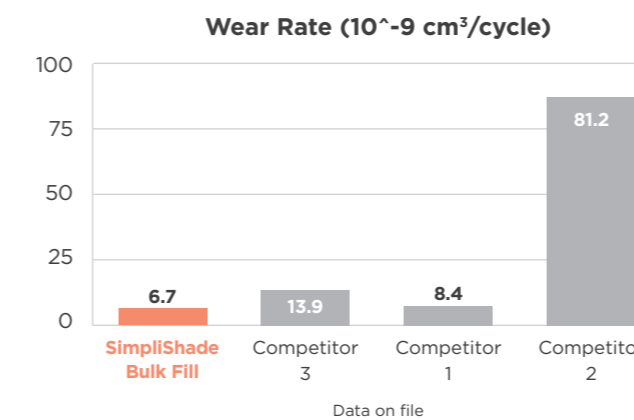

## Initial Gloss (GU)

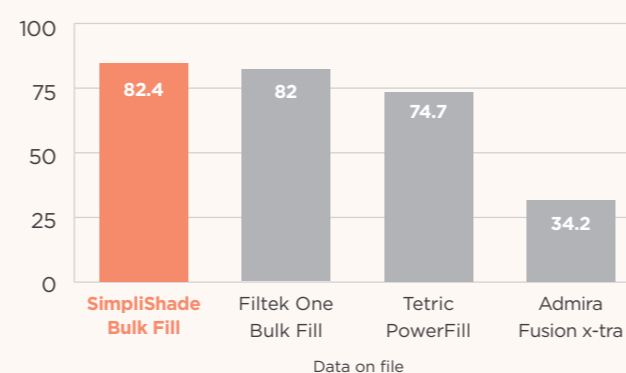


## Initial Gloss

SimpliShade Bulk Fill has comparable initial gloss to Competitor 3, and higher initial gloss than Competitor 1 and Competitor 2, resulting in a more esthetic restoration.




## Handling Properties

Adaptive rheological response technology reduces viscosity during manipulation for easier sculpting and anatomy creation.

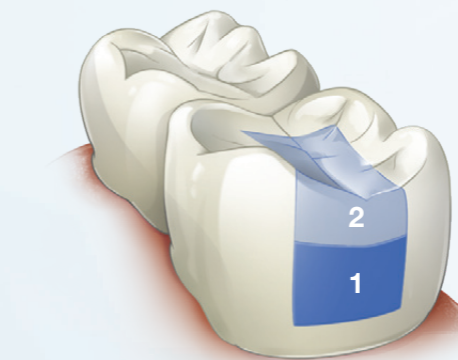

### Using flowable composite (1) and top up with SimpliShade Bulk Fill (2)

Enhance sealing and adaptation by applying flowable composite to line cavity (up to 2 mm.) Then simply and effortlessly fill cavities with SimpliShade BF up to 5 mm.

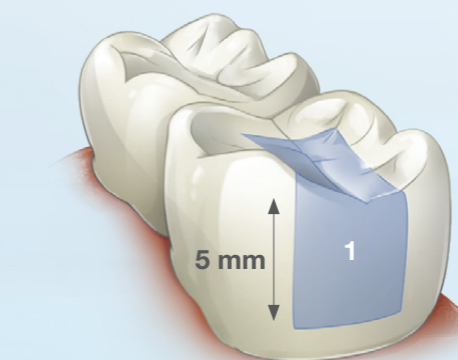

### Using only Bulk Fill for 5 mm

Simply and effortlessly fill cavities with SimpliShade BF up to 5 mm.

